

AM7515MT8A



Marca: Dino-Lite

Código del producto: AM7515MT8A

Short Description

5 Megapixel resolution

700-900x Magnification

Coaxial and brightfield illumination

Flexible LED Control (FLC)

Automatic Magnification Reading (AMR)

Descripción

The AM7515MT8A is a 5MP Dino-Lite Edge Series microscope that has a magnification range of 700-900x with built-in coaxial illumination. By using the Flexible LED Control (FLC), this model is capable of switching between or mixing, brightfield and coaxial illumination.



USB 2.0



Measurement
functionality



Adjustable
~700-900x



Flexible
LED control



Automatic
Magnification
Reading (AMR)



Magnification
lock




5 megapixel
2592 x 1944

Metal
housing



8 White LEDs
switchable

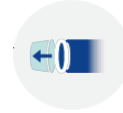

Exchangeable
cap

The coaxial illumination technique reveals details that are very difficult to see under normal light, for example when inspecting wafer plates, microchips or other microelectronics. The AM7515MT8A includes the Automatic Magnification Reading (AMR) feature which automatically detects and displays the magnification in the included Dino-Lite software. These unique features make the Dino-Lite AM7515MT8A a great inspection tool for material analysis, electronics inspection, or any similar application that requires high magnification, coaxial illumination, versatility and mobility.

The Dino-Lite AM7515MT8A is bundled with the user-friendly DinoCapture 2.0 software for Windows. For this model it includes functions such as Automatic Magnification Reading (AMR), Flexible LED Control (FLC), calibration, measurement, capturing & annotating images, and recording video. When focusing at such high magnification it is recommended to use a high-precision stand. The Dino-Lite RK-10A for example is a great add-on, it is a sturdy and stable high-end stand solution constructed of resilient stainless steel and lightweight aluminum and offers a very precise fine-focus adjustment.

Working distance/field of view/depth of field

MAGNIFICATION RATE	WORKING DISTANCE*	FIELD OF VIEW(X)	FIELD OF VIEW(Y)	DEPTH OF FIELD
700	6.4	0.56	0.400	0.01
750	6.3	0.52	0.375	0.01
800	6.1	0.48	0.350	0.009
850	6.0	0.46	0.330	0.009
900	5.9	0.43	0.315	0.008
Listed values may differ slightly	*Without front cap			Unit = mm

Especificación

Lighting	
Light/ LED type	White
Number of LEDs	8
LED on/off switchable:	Yes
Infrared filter	IR cut-filter >650 nm
Diffuser available	No
Emission filter	No
Polarizer	No
Optics	
Magnification	700-900x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	5 Megapixel (2592x1944)
Maximum frame rate	30 fps (max 20fps video recording)
Compatibility	
Interface	USB 2.0
Operating system	Windows 7, 8, 10 & 11, MacOS 10.12 and up
Included software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats (MacOS)	MOV (max 1.3MP)
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Metal housing
Magnification lock	Yes
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	140g
Cable length	1.8m
Features	
Special feature	Automatic Magnification Reading (AMR),

	Flexible LED Control (FLC), Coaxial illumination
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	Yes
Information	
Package contents	Microscope, Carry pouch, software CD, user manual, quick guide, calibration target, front cover box
Warranty information	2 years European warranty
Regulatory approval	CMOS, UKCA – available upon request – marketing@dino-lite.eu
Price range	€1000,00 - €1300,00

Product Gallery