

AM4115T-YFGW

Marca: Dino-Lite

Código del producto: AM4115T-YFGW



Short Description

1.3 Megapixel resolution
525nm excitation
570nm emission filter
20-220x magnification
To observe Cy3, TRITC

Descripción

The Dino-Lite AM4115T-YFGW digital microscope is optimized for research and viewing fluorescent objects by using 525nm LEDs. It has a 570nm emission filter that is designed to observe but not limited to Cy3 and TRITC.



USB 2.0



Measurement
functionality



Adjustable
~20x-220x



Magnification
lock



1.3 megapixel
1280 x 1024



Fluorescent
LEDs ~ 525nm



Standard working
distance

Compared to the traditional fluorescence microscope's band-pass type of emission filters, the Dino-Lite's long-pass type emission filter provides visibility and sensitivity over a


larger range of the fluorescence wavelength. The AM4115T-YFGW has the capability of switching the light source from the 525nm green excitation LED's to the white LED, which is convenient for locating the object and obtaining an easy focus.

» *Special models available on request.*

Working distance/field of view/depth of field

Magnification Rate Listed values may differ slightly	Working Distance * * Without front cap	Field of View(x)	Field of View(y)	Depth of Field Unit = mm
20	60.2	19.5	14.6	2.5
30	33.5	13.0	9.7	1.8
40	20.9	9.8	7.3	1.5
50	13.9	7.8	5.8	-
60	9.7	6.5	4.8	-
70	7.1	5.6	4.2	1.0
80	5.5	4.9	3.6	-
90	4.5	4.3	3.2	-
100	4.1	3.9	2.9	-
110	4.0	3.6	2.7	-
120	4.1	3.3	2.4	-
130	4.5	3.0	2.2	-
140	5.0	2.8	2.1	-
150	5.6	2.6	1.9	-
160	6.3	2.4	1.8	-
170	7.1	2.3	1.7	-
180	8.0	2.2	1.6	-
190	8.9	2.1	1.5	-
200	9.9	2.0	1.5	-
210	10.9	1.9	1.4	-
220	11.9	1.8	1.3	0.1

Special Light

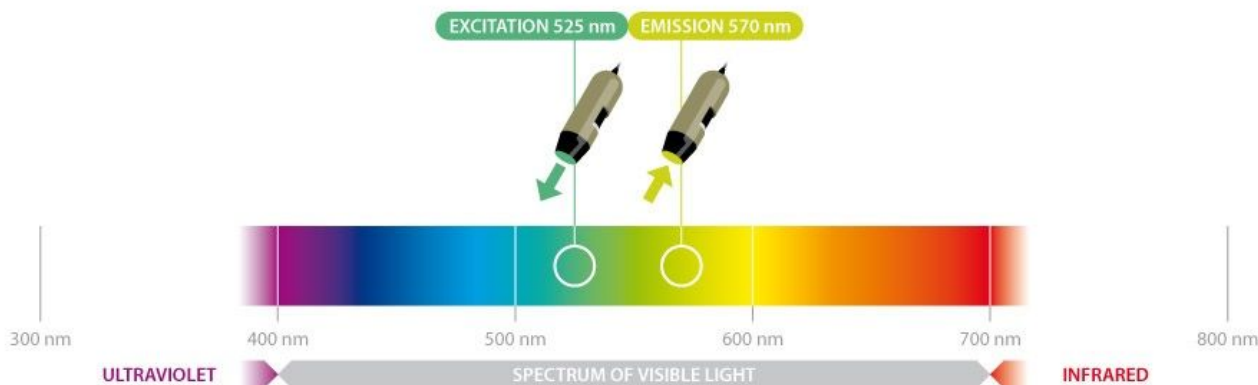

Especificación

Lighting	
Light/ LED type	525nm excitation + white
Number of LEDs	7 FL + 1 White
LED on/off switchable:	Yes
Infrared filter	No
Diffuser available	No
Emission filter	570nm, LP
Polarizer	No
Fluorophores	Cy3, TRITC
Optics	
Magnification	20-220x
Macro zoom	No
Working distance	Standard
Lens type	Glass with anti-reflection coating
Sensor	
Sensor type	CMOS
Resolution	1.3 Megapixel (1280x1024)
Maximum frame rate	30 fps
Compatibility	
Interface	USB 2.0
Operating system	Windows 7, 8, 10 & 11, MacOS 10.9 and up
Included software	DinoCapture 2.0 (Windows), DinoXcope (Mac OS)
Supported image formats (Windows)	BMP, GIF, PNG, JPG, TIF, RAS, PNM, TGA, PCX, MNG, WBMP, JP2, JPC, PGX
Supported video formats (Windows)	WMV, FLV, SWF
Supported image formats (MacOS)	JPEG, PNG
Supported video formats	MOV

(MacOS)	
Imaging standards	DirectShow, UVC
Wifi	Wireless-ready, requires the WF-10 WiFi streamer (optional)
Housing	
Housing material	Composite/ plastic housing
Magnification lock	Yes
Dimensions	10.5cm (L) x 3.2cm (D)
Weight	105g
Cable length	1.8m
Features	
Special feature	No
Measurement	Yes
Calibration	Yes
Microtouch sensor	Yes
ESD safe	No
Information	
Package contents	Microscope, carry pouch, software CD, calibration target, user manual
Warranty information	2 years European warranty
Regulatory approval	CE, FCC, ROHS
Price range	€700,00 - €900,00

Product Gallery

